BATAVIA

# MAXXGUN PRO

**Cordless massage tool**

Akku Massage-Gerät

Outil de massage sans fil

Accu massage tool

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Model: BT-MMT001  
Item-No.: 7063694  
[www.batavia.eu](http://www.batavia.eu)

### Overview

1. Arm rotation button
2. Battery power indicator light
3. Charging port
4. Gear indicator
5. Speed selection button
6. On/Off switch

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---

## Contents

1. Explanation of the symbols . . . . .	3
2. Massage tool specific safety warnings . . . . .	3
3. Safety instructions for battery powered devices . . . . .	3
4. Additional safety instructions for batteries and chargers . . . . .	4
5. Before the first use . . . . .	4
6. Contents of packing . . . . .	5
7. Intended use . . . . .	5
8. Instructions for use . . . . .	5
9. Technical data. . . . .	6
10. Cleaning and maintenance . . . . .	6
11. Disposal . . . . .	6
12. EC-Declaration of conformity . . . . .	6

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## Dear customer

Please familiarize yourself with the proper usage of the device by reading and following each chapter of this manual, in the order presented. Keep these operating instructions for further reference.

### Please read all safety instructions!

This manual contains important information on operating and handling the device. Keep these operating instructions for further reference. Please pass it on along with the unit if it is handed over to a third party.

## 1. Explanation of the symbols

The following symbols are used in the user manual or on the product:

**Warning!** To reduce the risk of injury, the user must read the instruction manual.

**Electric power tool, protection class II.** This sign symbolises that this unit corresponds to Protection

**Class II.** This means that the unit is equipped with either reinforced or double insulation between the mains circuit and output voltage or metal housing respectively.

**CE** stands for "Conformité Européenne", which means "In accordance with EU Regulations".

With the CE marking, the manufacturer confirms that this Electric tool complies with the applicable European directives.

**Do not dispose of electric power tools with domestic refuse.**

## 2. Massage tool specific safety warnings

- To reduce the risk of electric shock, fire and personal injury, this product must be used in accordance with the following instructions.
- Only allowed for adults.
- Apply pressure that is comfortable for you. Let the massager do the work. You should not be sore as a result of use.
- Use high speed and high pressure only for soft tissue parts of the body without causing pain and discomfort.

- Do not use high speed and high pressure on the head and any bone parts of the body.
- According to the massage needs of different parts please choose the right massage head.
- Frequent use of high-speed or high pressure to apply pressure on the same part will cause abrasions. If you feel pain and uncomfortable, please stop using the massage gun immediately.
- Do not get hair caught in the massage bit arm or any other moving part on the massager.
- Do not use massager on your head or near your genitals.
- Do not use on or near eyes or other highly sensitive areas.
- The equipment should not be put into the water.
- Do not discard and misuse this device.
- Carefully check your device before each use.
- Do not make any private modifications to the equipment.
- Do not run the device and charge the device in an unmanned state.
- Do not remove the screws by yourself.
- Turn off the device after one hour of operation and use it after rest for half an hour.
- The use of this device for an injured part of the body is prohibited without consulting the doctor.
- Do not block the vents of the motor.
- Never leave the appliance unattended as long as it is in operation.

**ATTENTION:** Do not operate the appliance at temperatures below 0 °C.

## 3. Safety instructions for battery powered devices

- 3.1 Remove the charger from the device after charging and do not use the device during charging.
- 3.2 Do not incinerate a discarded battery. It can explode in a fire.

- 3.3 Do not use an extension cord. Plug the charger cord directly into an electrical outlet.
- 3.4 Charge the battery only with the charger supplied with this cordless massager.
- 3.5 Do not use the charger in wet or damp conditions. It is intended for indoor use only. Do not use the charger near sinks or bathtubs. Do not immerse the charger in water.
- 3.6 Do not allow the cord to hang over the edge of a table or counter or touch hot surfaces.
- 3.7 The charger should be placed away from sinks and hot surfaces.
- 3.8 Do not use charger to charge any batteries other than the cordless massager battery. Other batteries may explode.
- 3.9 Do not operate charger if cord or plug is damaged. Replace damaged cord and plug immediately.
- 3.10 Do not operate charger if it has received a sharp impact, been dropped or otherwise damaged in any way. Have a qualified technician examine the charger and repair it if necessary. Do not disassemble the charger.
- 3.11 For best results, charge the battery when the temperature range of the battery and the work area is 0–40 °C (32–104 °F).
- 3.12 Unplug the charger when not in use and before cleaning or maintenance.
- 3.13 To reduce the risk of damage to the electric cord or plug, never carry the charger by the cord or yank the cord to unplug.
- 3.14 Always grasp the plug and pull to disconnect.
- 3.15 Always keep the cord away from heat, oil and sharp edges.

## 4. Additional safety instructions for batteries and chargers

### 4.1 Batteries

- Never attempt to open for any reason.
- Do not expose the battery to water.
- Do not expose the battery to heat.

- Do not store in locations where the temperature may exceed 40 °C.
- Keep batteries out of the reach of children.
- Charge only using the charger provided with the appliance/tool.
- When disposing of batteries, follow the instructions given in the section "Disposal".
- Do not charge damaged batteries.
- Under extreme conditions, battery leakage may occur. When you notice liquid on the batteries, proceed as follows: Carefully wipe the liquid off using a cloth. Avoid skin contact. In case of skin or eye contact, follow the instructions below.

**WARNING:** Never short-circuit the battery or expose it to moisture. Do not store it together with any metal parts which could short-circuit the contacts. The battery could overheat, start burning or explode.

**WARNING:** The battery fluid may cause personal injury or damage to property. In case of skin contact, immediately rinse with water. If redness, pain or irritation occurs seek medical attention. In case of eye contact, rinse immediately with clean water and seek medical attention.

### 4.2 Chargers

Your charger has been designed for a specific voltage. Always check that the mains voltage corresponds to the voltage on the rating plate.

#### **WARNING:**

- When not in use and before cleaning and servicing, disconnect the battery charger from the power supply.
- Never attempt to replace the charger unit with a regular mains plug.
- Never attempt to charge non-rechargeable batteries.
- The appliance/tool/battery must be placed in a well ventilated area when charging.
- The charger is intended for indoor use only.

## 5. Before the first use

Remove the machine and all accessories from the package. Keep the packaging

materials out of the reach of small children. There is a danger of suffocation!

## 6. Contents of packing

- 1× Maxxgun
- 1× Round head
- 1× Spinal head
- 1× Flat head
- 1× Precision point head
- 1× Nylon carry bag
- 1× Battery
- 1× Charger
- 1× Instruction manual

## 7. Intended use

This massage gun is intended for massaging soft body parts. This battery powered massage gun should be used for this only.

## 8. Instructions for use


### 8.1 Use

- a. To turn the device on, use the on/off switch (6) on the back of the device (See Fig. B).
- b. Put it in the desired speed using the speed selection button (5) above the on/off switch (6) (See Fig. B).
- c. Contact your body with the massager and let it float across your body. Get used to the feeling of vibration therapy.
- d. After the massage is complete turn the on/off switch (6) to the "OFF position and the light will turn off.

**NOTE:** The best process is to start the massager away from your body, then once it is operating, contact your body with the massager bit and allow the massager to float across your tissue. Once you have a feel for the pressure you desire, you can increase the pressure slowly until you reach a comfortable balance.

### 8.2 Battery indicator

Once the massager is switched on, the battery power indicator light (2) will indicate the battery charge level. When there is only one light on, recharge your massager (See Fig. A).

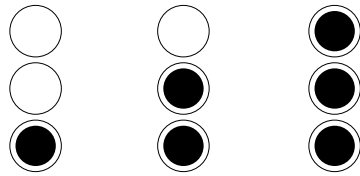Recharge    Partially charged    Fully charged

### 8.3 Adjusting the speed

To adjust the speed of the massage head, use the speed selection button (5) on the back of the device. Click it continuously to cycle through the speed gears. Left light gives 1800 rpm, middle light 2100 rpm and right light 2400 rpm.

### 8.4 Adjusting the arm

- a. Push and hold the arm rotation button (1) (See Fig. A).
- b. While holding the button, rotate the arm to 1 of 3 positions.
- c. Once you hear the arm lock in place, you can release the arm rotation button (1).

**NOTE:** According to the different needs of the body, select the appropriate massage head to insert into the round hole in front of the massage gun.

### 8.5 Charging

- a. Insert the charger plug into the charging port (3) on the bottom of the massager (See Fig. A).
- b. Plug the charger into a wall receptacle. When fully inserted, the three battery indicator lights on the massager handle will flash one by one alternatively, and the battery will start to be charged.
- c. When one of the three battery indicator lights (2) turns into solid bright while the other two continue flashing one by one alternatively,  $\frac{1}{3}$  of the battery is fully charged (See Fig. A).
- d. When two of the three lights turn into solid bright, meanwhile the third light continues flashing,  $\frac{2}{3}$  of the battery are fully charged. Once the three lights turn into solid bright, the battery is fully charged.

## 9. Technical data

Battery:	24 Volt Li-Ion
Motor power:	250 W
Stroke:	12 mm
Stall force:	30 kg
Speed:	1800/2100/2400 rpm
Frequency:	
1st gear	30 Hz
2nd gear	35 Hz
3rd gear	40 Hz
Weight:	1 kg
Charging time:	2–3 hrs
Noise:	60 dbA

## 10. Cleaning and maintenance

- Be sure to unplug the unit and allow it to cool before cleaning.
- Wipe the surface of the device with a slightly damp towel and dry with a soft cloth.
- Store the machine and the accessories in a dry place.

## 11. Disposal

*Do not dispose of electric power tools with domestic refuse.*

The electric power tool is shipped in packaging to reduce transport damage.

This packaging is a raw material and as such can either be reused or can be fed back into the raw material cycle. The electric power tool and its accessories are made from various materials such as metals and plastics. Take defective components to a special refuse collection point. Ask about these at your specialist shop or local council.

### 11.1 Batteries

Think of the environment when disposing of batteries. Contact your local authorities to find out where your nearest disposal area is. Do not dispose of batteries with domestic refuse. Do not dispose of by burning, risk of explosion.

## 12. EC-Declaration of conformity

We, the **Batavia B.V., Weth. Wassebaliestraat 6d, NL-7951 SN Staphorst**, declare by our own responsibility that the product **Cordless massage tool, Item-No 7063694, Model Nr. BT–MMT001** is according to the basic requirements, which are defined in the European Directives **Electromagnetic Compatibility 2014/30/EU (EMC), RoHS Directive 2011/65, Appendix II Amendment Directive (EU) 2015/863** and their amendments. For the evaluation of conformity, the following harmonized standards were consulted:

**EN55014-1:2017**

**EN55014-2:2015**

Staphorst, 1 december 2019

Meino Seinen, QA Representative  
Batavia B.V., Weth. Wassebaliestraat 6d,  
7951 SN Staphorst, Netherlands

The product and the user manual may be subject to changes. Technical data may be changed without prior notice.

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# BATAVIA

## 2 YEAR WARRANTY

### This product has got a 2 year warranty

Dear Client, if for any reason this product is not working, please ensure you contact our Client Service Centre. Ensure you have your original receipt of purchase. This warranty covers all defects in workmanship or materials in this Batavia product for a two year period from the date of purchase. The warranty does not cover any malfunction, or defect resulting from misuse, neglect, alteration, or repair.

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Meino Seinen, Batavia B.V.,